

Otto Bock[®]

QUALITY FOR LIFE



Start^{M1} Intro

Product Information





Start M1 Intro



The Start M1 Intro combines the concepts of modularity, quality and availability like no other product in its class. It assumes the place of a readily available fleet chair within the Start family, leaving more complex configurations to the specialised models in the product line. 20 versions of the Start M1 Intro are available from inventory in order to readily meet short-term demand – even for large numbers of units. Key elements like the segmented swing-away footrest, removable side panels and drive wheels with quick-release axles are included as standard equipment. The seat depth and back height can be adjusted to the required dimensions at any time. Selected accessories are included with the shipment for easy installation.

Features

- Round aluminium tube frame
- Tip-assist integrated into the frame
- Segmented footrest, angle-adjustable
- Swing-away removable footrest holder
- Height-adjustable side panel, with arm rest
- Nylon seat and back upholstery, black
- 8" PU steering casters
- 24" drive wheels with quick-release axles
- Knee lever wheel lock
- Frame colour silver metallic

Accessories

- Transport wheels
- Anti-tipper (single or double)
- Crutch holder
- Various side panels
- Tray
- Elevating leg rests
- Lap belt
- Seat cushion
- Spoke protectors
- Motor vehicle transportation kit



Steering caster 8" x 50 mm PU



Tip-assist



Knee lever wheel lock



Drive wheel with aluminium push ring




Version with drum wheel lock



Version with back angle adjustment

Versions

Seat width	Drive wheel	Back	Wheel lock	Factory setting
40.5 cm	24" x 1 3/8" PU	Standard 42.5 cm	Knee lever wheel lock	Seat width 43 cm, seat height 49 cm, 3°
43 cm	24" x 1 3/8" PU	Standard 42.5 cm	Knee lever wheel lock	Seat width 43 cm, seat height 49 cm, 3°
45.5 cm	24" x 1 3/8" PU	Standard 42.5 cm	Knee lever wheel lock	Seat width 43 cm, seat height 49 cm, 3°
48 cm	24" x 1 3/8" PU	Standard 42.5 cm	Knee lever wheel lock	Seat width 43 cm, seat height 49 cm, 3°
50.5 cm	24" x 1 3/8" PU	Standard 42.5 cm	Knee lever wheel lock	Seat width 43 cm, seat height 49 cm, 3°
40.5 cm	24" x 1 3/8" PU	Standard 42.5 cm	Drum wheel lock	Seat width 43 cm, seat height 49 cm, 3°
43 cm	24" x 1 3/8" PU	Standard 42.5 cm	Drum wheel lock	Seat width 43 cm, seat height 49 cm, 3°
45.5 cm	24" x 1 3/8" PU	Standard 42.5 cm	Drum wheel lock	Seat width 43 cm, seat height 49 cm, 3°
48 cm	24" x 1 3/8" PU	Standard 42.5 cm	Drum wheel lock	Seat width 43 cm, seat height 49 cm, 3°
50.5 cm	24" x 1 3/8" PU	Standard 42.5 cm	Drum wheel lock	Seat width 43 cm, seat height 49 cm, 3°
40.5 cm	24" x 1 3/8" PU	Angle-adjustable, 47.5 cm	Knee lever wheel lock	Seat width 43 cm, seat height 49 cm, 3°
43 cm	24" x 1 3/8" PU	Angle-adjustable, 47.5 cm	Knee lever wheel lock	Seat width 43 cm, seat height 49 cm, 3°
45.5 cm	24" x 1 3/8" PU	Angle-adjustable, 47.5 cm	Knee lever wheel lock	Seat width 43 cm, seat height 49 cm, 3°
48 cm	24" x 1 3/8" PU	Angle-adjustable, 47.5 cm	Knee lever wheel lock	Seat width 43 cm, seat height 49 cm, 3°
50.5 cm	24" x 1 3/8" PU	Angle-adjustable, 47.5 cm	Knee lever wheel lock	Seat width 43 cm, seat height 49 cm, 3°
40.5 cm	24" x 1 3/8" PU	Angle-adjustable, 47.5 cm	Drum wheel lock	Seat width 43 cm, seat height 49 cm, 3°
43 cm	24" x 1 3/8" PU	Angle-adjustable, 47.5 cm	Drum wheel lock	Seat width 43 cm, seat height 49 cm, 3°
45.5 cm	24" x 1 3/8" PU	Angle-adjustable, 47.5 cm	Drum wheel lock	Seat width 43 cm, seat height 49 cm, 3°
48 cm	24" x 1 3/8" PU	Angle-adjustable, 47.5 cm	Drum wheel lock	Seat width 43 cm, seat height 49 cm, 3°
50.5 cm	24" x 1 3/8" PU	Angle-adjustable, 47.5 cm	Drum wheel lock	Seat width 43 cm, seat height 49 cm, 3°

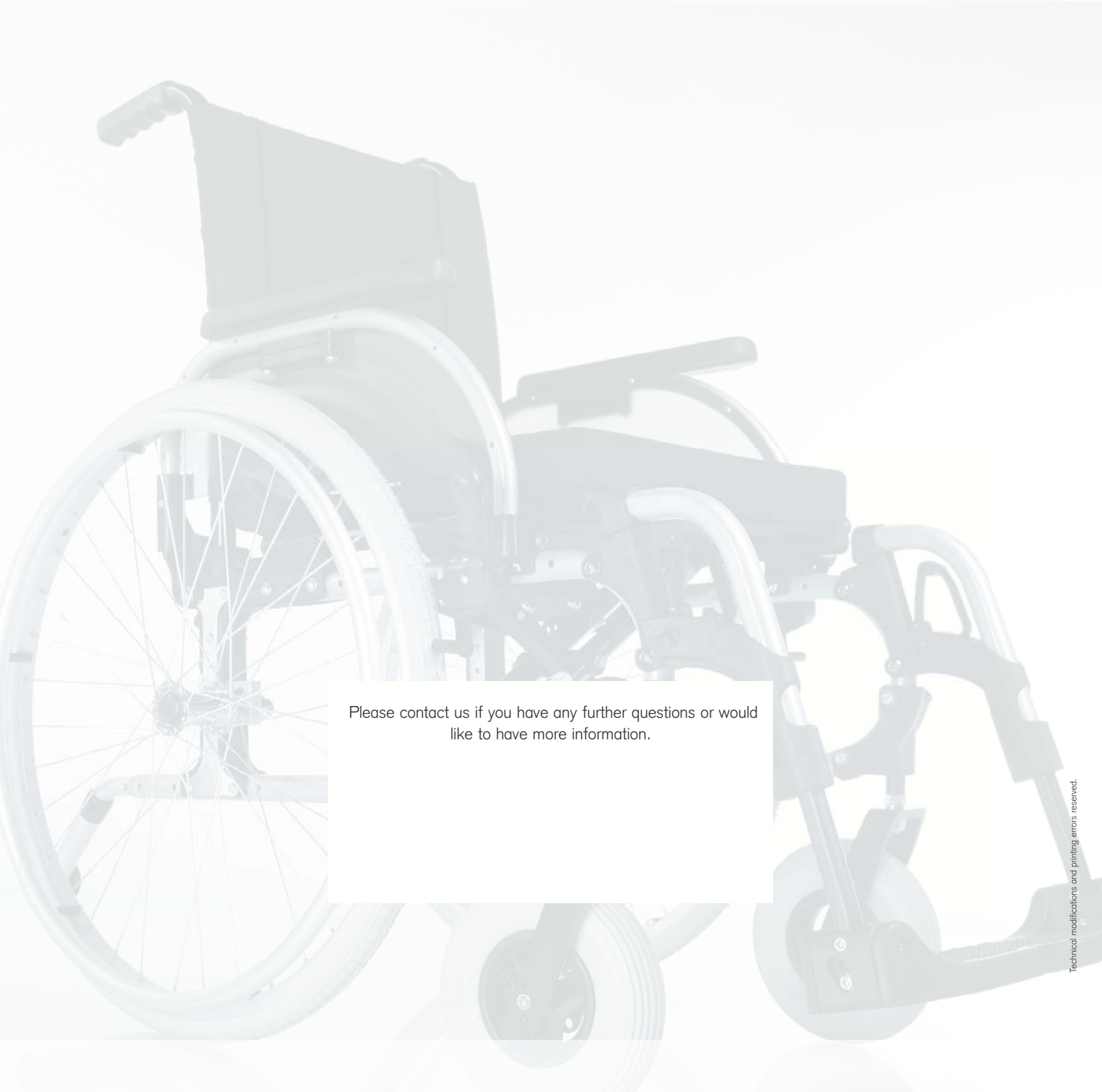

Please contact us if you have any further questions or would like to have more information.

Technical modifications and printing errors reserved.

Otto Bock®

QUALITY FOR LIFE

Otto Bock HealthCare GmbH

Max-Näder-Straße 15 · 37115 Duderstadt/Germany · Telephone +49 5527 848-1590 · Telefax +49 5527 848-1676 · reha-export@ottobock.de · www.ottobock.com

© Otto Bock · 646D36e-GB · 09.2009/1

